
Digital Personal Scale

USER'S MANUAL

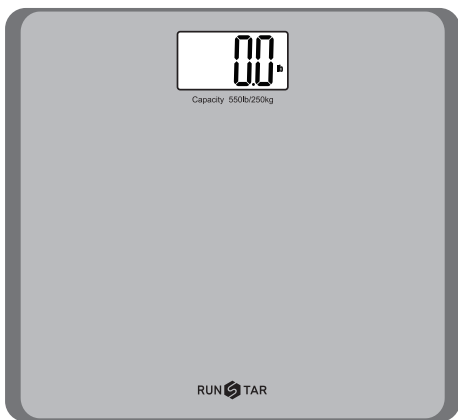Model: BG1911L

Contents

• ENGLISH	01-03
-----------	-------

• DEUTSCH	04-06
-----------	-------

• FRANÇAIS	07-09
------------	-------

• ITALIANO	10-12
------------	-------

• ESPAÑOL	13-15
-----------	-------


Thank you for choosing our electronic bathroom scale series.
To ensure proper function and to maximise the lifespan of your new scale, please read this manually carefully before use.

Key features

1. High precision strain gauge sensors system.
2. Tempered glass platform.
3. Capacity: 550lb/250kg
4. Division: 0.2lb/100g
5. LCD display
6. Auto on
7. Auto zero & Auto off
8. Low battery and over load indication
9. Free switch within three different units: kg/lb/st
10. Battery function: 3*AAA(1.5V) Batteries(included)

Operation

1. Remove the plastic strip if batteries are provided.

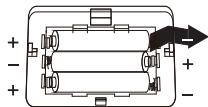2. Place the scale on a hard and flat surface.
A clean and hard floor is the best solution.

3. Step on the scale directly and
the weight will display.


4. Stay still on the platform. During weighting, indication will be flashing on the display. Wait till the indications stop flashing.
Your weight will be displayed for about 10 seconds, and then the unit will switch off automatically.
5. To switch units (kg-lb-st:lb): Press the "UNIT" button on the back of the scale display.

Caution

1. To ensure accuracy, we suggest you take measurement at a fixed time every day (e.g. 7:00pm-9:00pm).
2. Please change new batteries when "LO" appears on the LCD which indicates that the battery is out of power.

3. "Err" showed on LCD indicates overload .

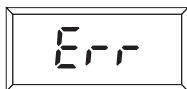Attention

1. Please take off shoes before taking measurement to avoid scraping or scratching the platform.
2. The platform may be pretty slippery when wet, thus please make sure both the platform and your feet are dry before using.
3. Please use soft tissue with alcohol or glass cleaner to clean the surface if it is dirty. No soap or other chemicals. Keep it away from water, heat and extreme coldness.
4. The scale is high precision measuring device. Never jump or stomp on the scale or disassemble it and please handle it carefully when moving it.
5. The scale is for family use only and not suitable for professional use .

FCC Statement

Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

FCC Statement NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help with the supplied fastener. DO NOT hang the unit by the power cord.

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