
SMART BODY FAT SCALE FG2016LB-B

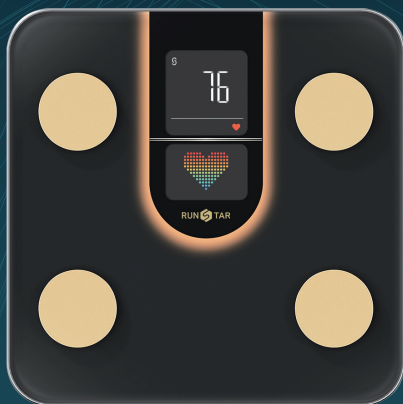User Manual

Contents

● ENGLISH	01-08
-----------	-------

● DEUTSCH	09-16
-----------	-------

● FRANÇAIS	17-24
------------	-------


● ITALIANO	25-32
------------	-------

● ESPAÑOL	33-40
-----------	-------

Bluetooth Scale Body Index

- Body weight
- BMI (Body mass index)
- Subcutaneous fat
- Body fat rate
- Muscle rate
- Fat-free body weight
- Heart rate
- Visceral fat
- Body water
- Skeletal muscle
- Muscle mass
- Bone mass
- Protein
- BMR
- Body age
- Cardiac Index
- Fat Mass
- Water Weight
- Protein Mass
- Ideal Body Weight
- Obesity Level
- Body Type

Trouble shooting

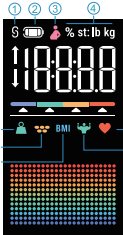1. Connected
2. Battery volume
3. Maternity mode
4. Units
5. Body Weight
6. Body fat rate
7. BMI
8. Heart Rate
9. Muscle rate

No display

Change battery

Change battery

Overload

Calibrating

For body weight: Place the scale on an even surface and wait at least 10 seconds before stepping on the scale. The reading should display 0.0 kg before stepping on.

During body fat testing process, the display shows bars in different colors in sequence, which refer to different body types.
After measured, the display shows the body fat rate and determinate color bar.

The correspondence between colors and body types:

light blue=thin


green=normal

orange=fat

red=obese

For example:

The below display shows the **body④ types** is fat.

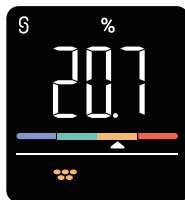
*Body type is determined based on BMI, with reference to the internationally uniform obesity classification standards recommended by the World Health Organization (WHO).

Caution

Conductive Areas

1. Do not shake or drop the scale.
2. Do not immerse the scale in water or use chemical cleaning agents. Clean the scale with a slightly damp cloth.

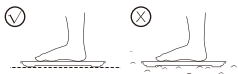3. Do not place the scale on any carpet or uneven flooring. Always place the scale on hard, flat, dry surface before measurement.

4. Always keep your scale in horizontal position.

5. Any measurement obtained using this device is for reference only and should not be considered as a medical opinion.

6. Consult your primary care doctor or physician before making changes to your diet, exercise plan or physical activity.

Specifications

Size: 10.94*10.94*1.06in/278*278*27mm

Weight: 3.17lb/1.44kg

Batteries: 4 x 1.5V AAA batteries

Capacity: 400 lb / 180 kg / 28 st

Division: 0.1lb/50g

Unit: lb/st/kg

Accuracy Range

50kg:±0.3kg	100kg:±0.4kg	180kg:±0.5kg
110lb:±0.66lb	220lb:±0.88lb	400lb:±1.1lb


Operating instructions

The Bluetooth Scale can connect with your mobile phone (Android 6.0+ & IOS 8.0+) or iPad via Bluetooth 4.0. Download the 'Starfit' app free from the App Store or Google Play.


Pairing the scale with your smartphone


Turn on Bluetooth on your smartphone.

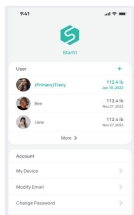Click [Account]→
[My Device] to enter
the device adding
page.

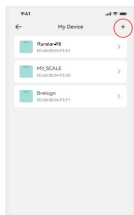
Power on the scale,
press[+] to select
"Runstar-R9" to
add the device.



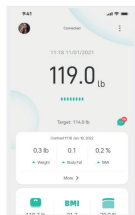
Come back to [Measure] page.
Step on the scale and keep standing
for 5-10 seconds with bare feet. The
weighing has finished after the
analysing data come out.

View reports, delete data, set reference data, share progress by recent, weekly, month and year.

Click [Chart] to view your measurement history. Click [User Data] to enter the calendar list and select a date to view or delete data. You can also compare data from different dates or share it with friends via Facebook, Instagram or SMS.



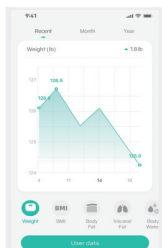
Click on the Body Composition icon to see a detailed description of the corresponding name and measurement.

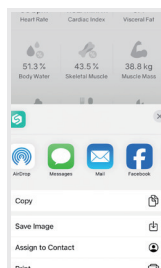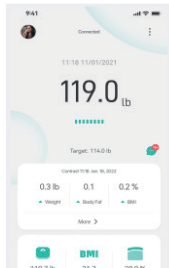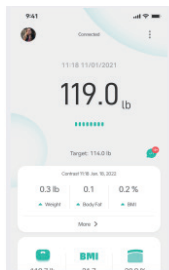Add/Delete Users & Basic Settings

Under [Account] you can add user(s) [+]. To delete an account, swipe the account to the left and select [Delete].

Set goals, Weight units, Weighing done sounds, Languages, Passwords can all be set under the [Account] page. Please note that you can change the weighing units by pressing the "Units" button (under the battery compartment cover), or you can also [Switch Units] from the [Account] page of App.

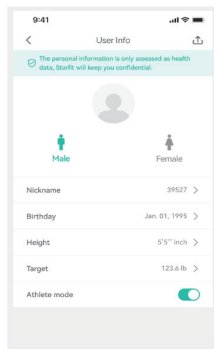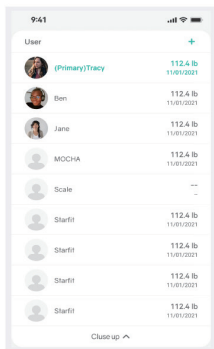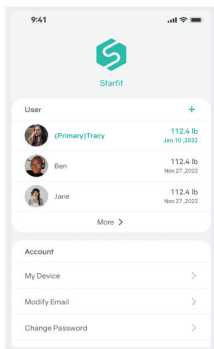Frequently Asked Questions&Answers

1.Why doesn't the scale work? Why do the data on screen disappear in a flash?

- Please check if the batteries are installed properly, replace batteries if necessary.

2. The scale cannot connect with App.

- Check if the phone's software is iOS 8.0 or higher or Android 6.0 or higher.
- Download and open the newest version of the Starfit App.
- Enable Bluetooth on your phone. For Android 6.0 or higher, location needs to be enabled.

3.No body fat measured when weighing.

- Step on with dry and bare feet.
- Ensure Bluetooth is enabled and working.
- Ensure personal data is entered.

Warning Note

Changes or modifications not expressly approved by the party responsible for compliance could avoid the user's authority to operate the equipment.

FCC Statement

Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference.
- (2) This device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC ID: 2AP3Q-FG2016LB-B

Model No.: FG2016LB-B

Guangdong Welland Technology Co., Ltd.

No.85, Minke East Road, Minying Science & Technology Park,
Shiqi District, Zhongshan, Guangdong, China

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Guangdong welland Technology Co. , Ltd.