Smart Full Body Composition Scale User Manual

FCC ID: 2AP3Q-FG2305ULB

Contents

● ENGLISH	1-10
-----------	------

● DEUTSCH	11-20
-----------	-------

● FRANÇAIS	21-30
------------	-------

● ITALIANO	31-40
------------	-------

● ESPAÑOL	41-50
-----------	-------

Visiting <https://vip.runstar.store> or scan the QR code below or to register and get lifetime warranty:

Thanks for choosing Runstar Smart body composition Scale. This scale can help you conveniently measure body weight, Partial body fat & Skeletal muscle, BMI, Body fat rate, Visceral fat, Body water, Skeletal muscle rate, Muscle mass, bone mass, Protein, BMR, Body age and etc.

IMPORTANT: Please read these instructions in this manual carefully before using the device. Failure to comply with the instructions given in the manual or using this scale in ways other than the ones mentioned in this manual may result in injury to yourself or damage to the scale.

Starfit APP Functions

- ◆ Weight
- ◆ Muscle weight
- ◆ Body age
- ◆ Protein rate
- ◆ Partial body fat
- ◆ Body fat
- ◆ Bone mass
- ◆ Muscle
- ◆ Subcutaneous fat
- ◆ Partial muscle rate
- ◆ Visceral fat
- ◆ Water weight
- ◆ Standard weight
- ◆ Height record
- ◆ Protein mass
- ◆ BMR
- ◆ Without fat
- ◆ BMI
- ◆ Degree of obesity
- ◆ BFR


Indication Symbols

LO	Low battery - Please recharge by Type C USB cable
Err	Overload (Scale will show Err,but app can not show this icon)
C	An error has occurred during measurement. Step off the scales and step onto it again to repeat the measuring process.

WHAT'S IN THE PACKAGE

- ◆ 1 X Smart body fat scale
- ◆ 1 X User manual
- ◆ 1 X USB C Cable

TFT Screen

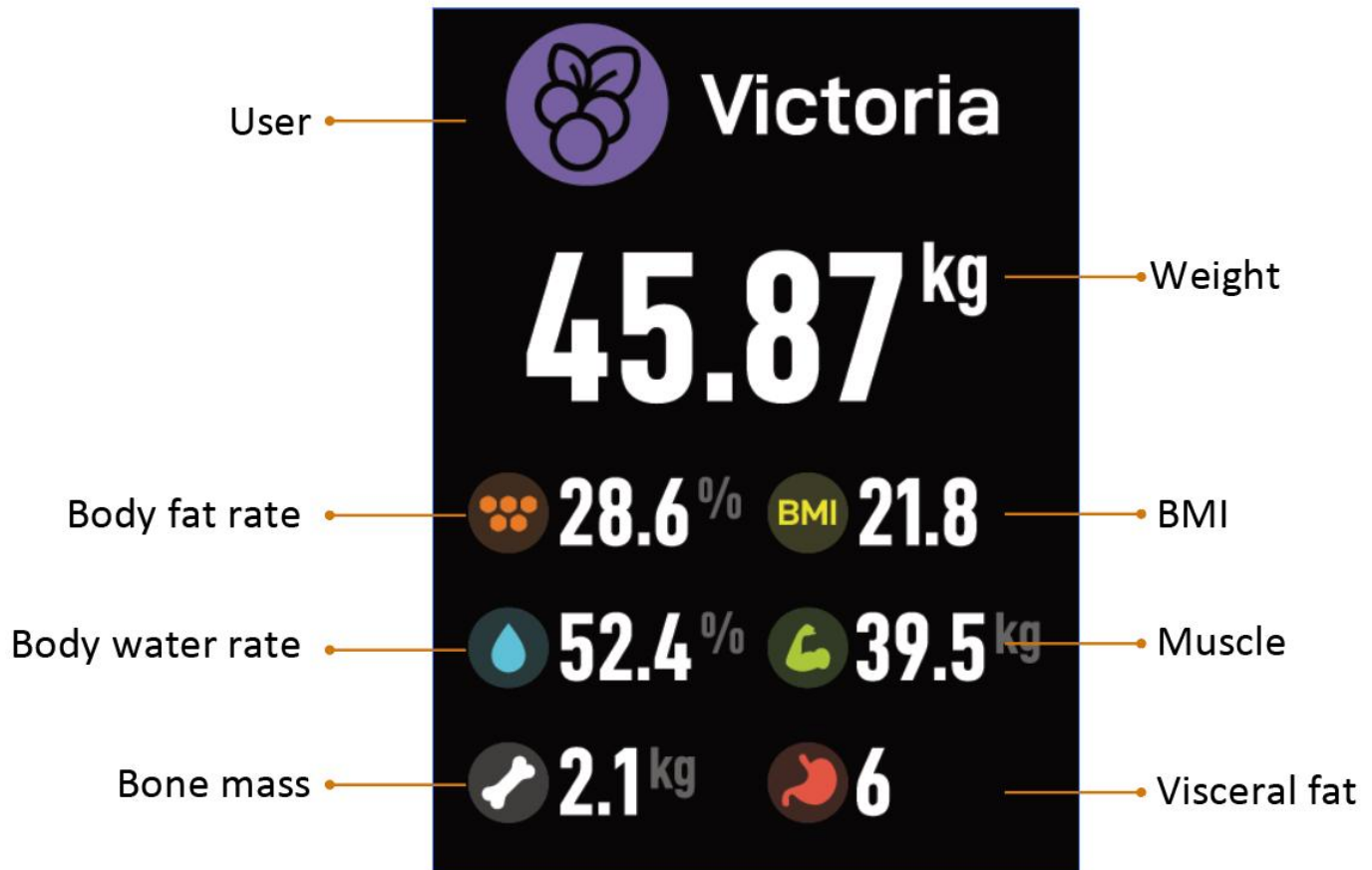WARNING

- ★ DO NOT use with medical implants such as pacemakers.
- ★ DO NOT stand on the edge of the scale or jump on it;
- ★ DO NOT overload the scale (Max. 396lb/180kg)
- ★ DO NOT drop scale or drop and objects on it as this may damage the sensors.
- ★ DO NOT immerse the scale in water or use chemical cleaning agents. Clean the scale with slightly damp cloth.
- ★ Not recommended for pregnant women.
- ★ Not recommended for infants, toddlers, and children under 10 years of age.
- ★ Any measurement obtained using this device is for reference only and should not be considered as a medical opinion.
- ★ Consult your primary care doctor or physician before making changes to your diet, exercise plan or physical activities.
- ★ Always place the scale on a hard, dry and flat surface before measurement.
- ★ Make sure your feet are dry before stepping on the scale.
- ★ Connect the batteries in accordance with the correct polarities indicated.
- ★ Remove batteries if the scale is not be use for a prolonged period of time.

SPECIFICATION

Weight unit: lb / kg / st:lb (switch the units via APP setting)

Weight limit: 11lb-396lb / 6-180kg

Power Supply: 500mAh rechargeable battery

Account number: No limitation (each account can add 24 persons)

Operating temperature: 10-40°C

Operation humidity: 20%-90%

Support Devices:

iOS 8.0 or higher

Android 6.0 or higher

Install Starfit App

1. Search “Starfit” in App Store or Google Play or scan the QR Code below:

2. Download and install the App on your device.
3. The App icon will appear on your phone after the installation is completed.


Special selling points:

1. This model features a colorful TFT display and offers offline functionality. Setting up the scale and connecting it with the app is crucial for the initial use to analyze body fat and other parameters. However, after the initial setup, you can conveniently measure body fat, BMI, and body types without the need to connect with the app every time. The TFT display on the scale allows you to view these measurements. It's important to note that offline measurement will not be activated if there is a weight difference greater than 2kg compared to the last measurement.
2. During offline mode, any successful measurement results will be automatically synchronized with your app during the next connection between the scale and the app. This ensures that your data is seamlessly updated and accessible for further analysis and tracking.

How to use Starfit APP


1. Press the button besides USB C charging port to turn the scale on.
2. Place the scale on the hard flat floor without carpet.
3. Turn on Bluetooth on your smart phone.

4. Register your own ID in Starfit App using email.

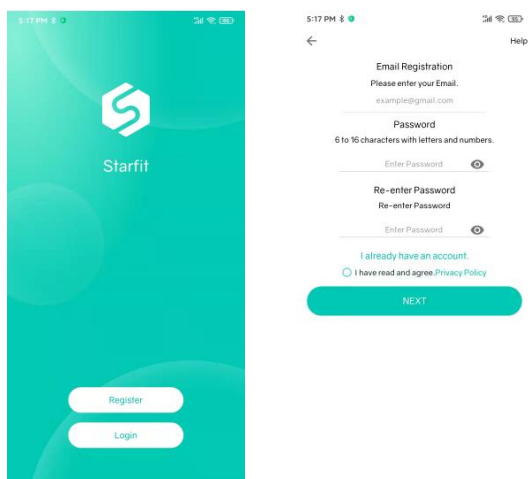
5. Add personal data and confirm OK;

Note: Entering incorrect gender, age, or height information may result in inaccurate measurement results

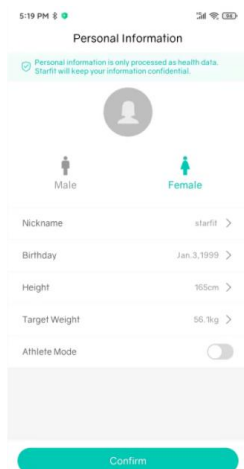
6. Pair the scale through Bluetooth. Click 'Account' → 'My Device' → '+' → 'Searching Bluetooth for pairing' to find the scale (Keep scale on) to pair. And now your phone is connected with the scale;

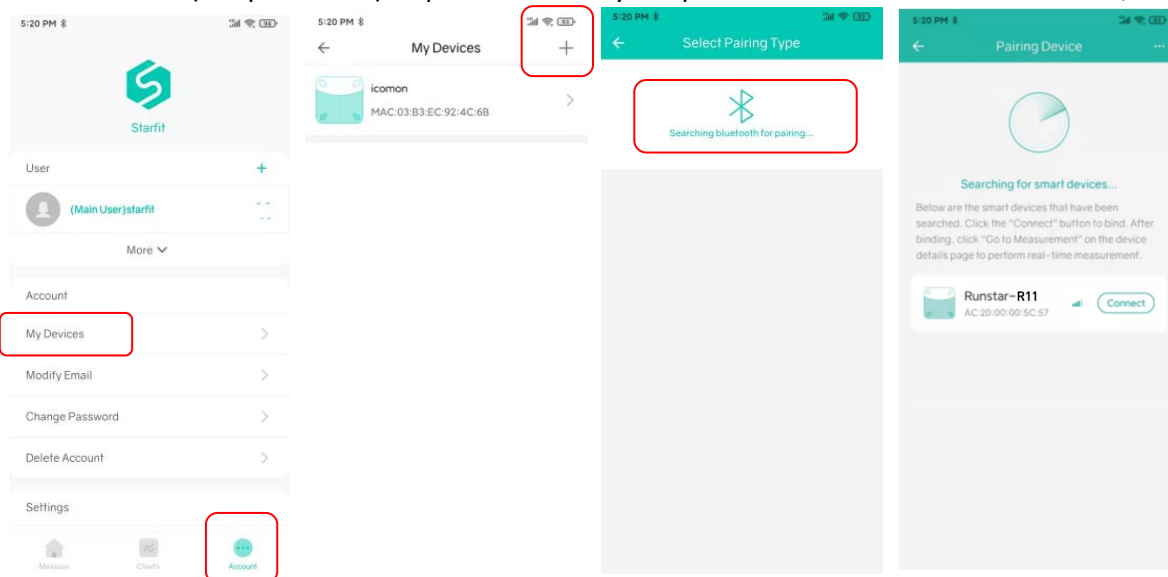
After connected, **S** will be displayed on the scale screen, and "Connected" will be displayed at the top of the App's "Measure" page. If it displays "Disconnected", it means that the pairing failed, please repeat the above steps.

7. Step on and off the scale once to activate the scale. Wait until the display shows "0.00" (lb/kg/st). Keep Starfit App on 'Maesure' page.

8. ① Step on scale with bare feet. Please make sure that the soles of your feet are in full contact with the electrodes on the scale (pay attention to Left & Right direction), as shown in the figure below.

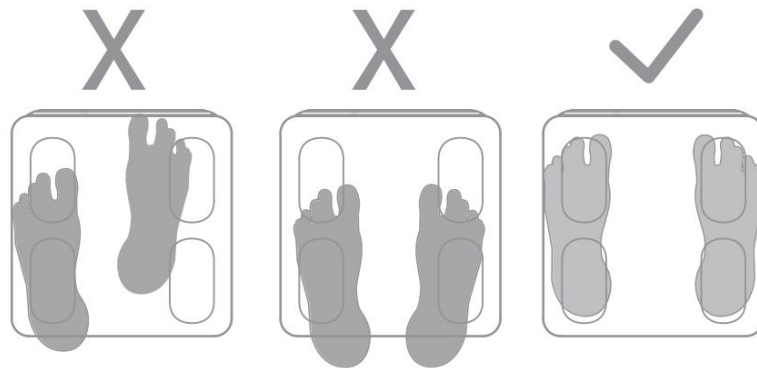② After standing firmly, grab the handle, making sure your thumbs touch the electrodes below the display screen, and the rest of your fingers fully touch the electrodes at the back of handle; stand up straight, spread your legs, stretch your arms naturally and keep them at a 45° angle to your body, as shown in the figure below.

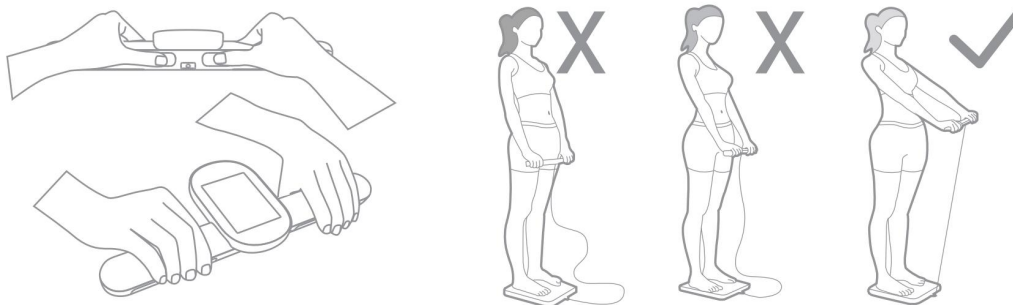Note:

If you are not barefoot or your feet are too dry, or your weighing posture is incorrect, only your weight and BMI can be measured.


9. The number on the display is flashing when the scale is weighing and analyzing.

Please stand still and keep same posture for about 20S;

App will show measure schedule and with music playing during body fat measurement.

After the measurement is completed, you can view all the results on the App.

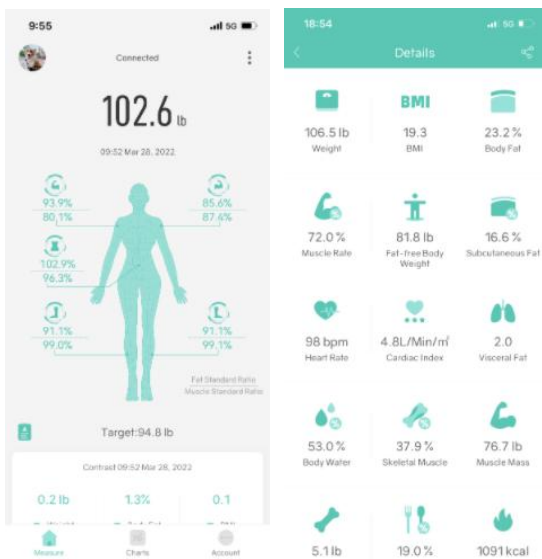At the same time, you can read 7 core data directly on TFT screen, including body fat rate, BMI, Weight, Body Water Rate, Bone Mass, Muscle, Visceral Fat.

Besides, after body data measuring, the screen will show the body data difference which compared with last body data measurement.

During body data measurement, scale display will show as below:

Body weighing

Measuring

Weight Comparison

Body fat rate


Fat Rate Comparison

Muscle


Muscle Comparison

Body Date Summary

10. You can click on 'More' to enter 'Details' to get the complete measurement results. At the bottom of the 'Details' page, you can perform 'Report printing'

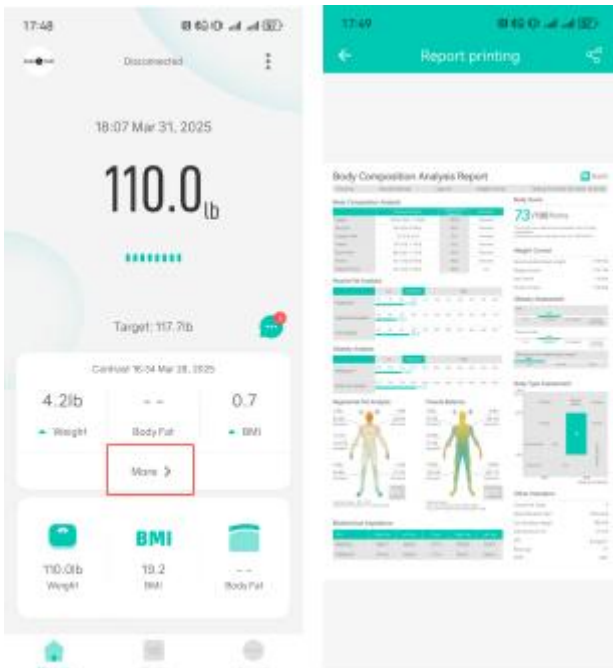
11. Click on 'Charts' to view past measurement records in the form of a curve chart. Click on 'User Data' to view past records of weight, BMI, BFR, muscle, bone mass, etc. You can also delete data or compare data from different dates.

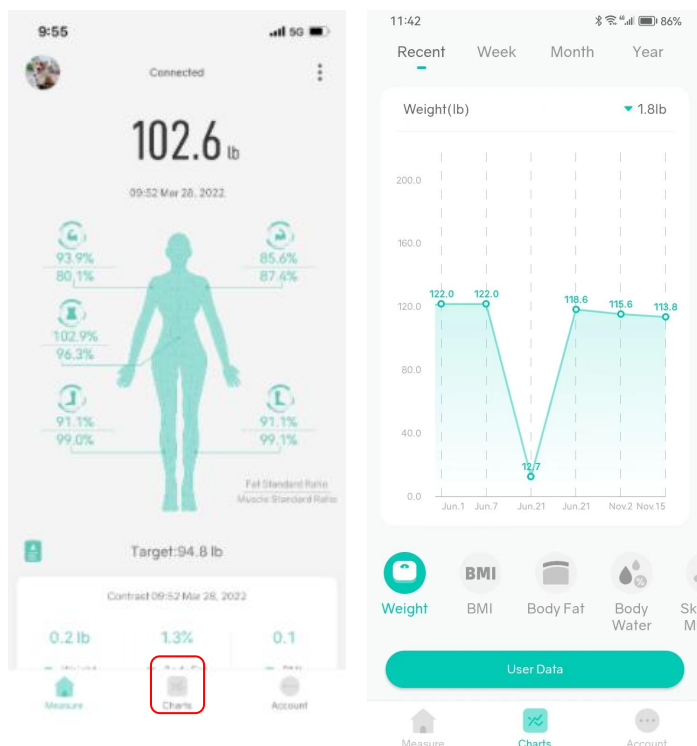12. Add or Switch user:

Add users by clicking the user avatar in the upper left corner of the 'Measure' page → 'Add User' or on the 'Account' page → '+'. Each account supports up to 24 users;

Switch user: Click the corresponding user avatar to switch users

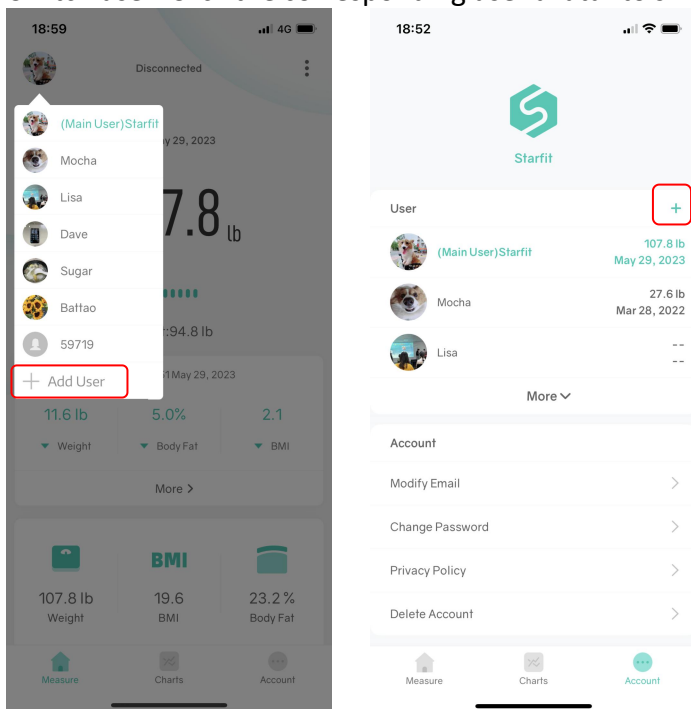13. When the battery is low power, TFT display will indicate and flick as below picture. Please charge the scale with USB C cable.

✂It is recommended to use a 5V/1A adapter.

TROUBLE SHOOTING GUIDE

1. Bluetooth not connected

- The scale and smart phone must be both have Bluetooth on(scale lights up means Bluetooth on).
- Check the iOS version of your smart phone, it must be iOS 8.0 or higher/Android 6.0 higher.
- Make sure your phone allows the APP to use Bluetooth. If not, please turn it on in the phone settings.

2. No body fat data measured when weighing

- User must be bare foot to measure the body fat.
- Do pay attention for correct handle bar grip posture during measurement.
- Make sure the soles of feet are slightly damp (apply some body lotion) and the scale surface is dry.
- The measurement process was not fully completed: Make sure to stand still on the scale for about 20s until the user's Avatar, Name and Weight appear on the display.
- The scale disconnects from the app.
- Scale is connected to another user's phone - scale shows connected, but the APP shows disconnected

3. The scale is producing inconsistent weight.

Check whether the scale is placed on a hard, flat floor without carpet.

Check each sensor foot on the back to make sure nothing is stuck to the bottom of it.

Press the 'reset' button and measure again.

4. The App is display the incorrect fat-muscle ratio

The body parameter setting is incorrect, check to make sure you entered the correct gender, height, and age.

5. How to recalibrate the scale?

Press the "Reset" button or step on the scale repeatedly so that it displays "0.00". If it still doesn't work, contact the customer service team for a replacement.

6. The scale does not turn on

Charge for 20mins and try again. If it doesn't work, contact the customer service team to get a replacement.

7. Scale only show weight & BMI?

- Check if Bluetooth is connected
- Check if the measurement posture is correct: whether the hands and feet are fully in contact with the 8 electrodes. This manual provides instructions for the correct measurement posture.
- Keep your legs and feet apart. The skin proximity of your legs and feet will affect the transmission of the measurement current.
- Stand still for about 20-30 seconds to complete the entire measurement.

FCC Statement

Operation is subject to the following two conditions:

(1) This device may not cause harmful interference.

(2) This device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.

- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Battery Advice

* Place the scale far away from the children in case of falling down or crashing.

Disposal

The packaging material is recyclable. Please dispose of it in an environmentally friendly manner and take it to a recycling center.

Dispose of the device in an environmentally friendly manner. Dispose of it at a recycling center for used electrical and electronic devices.

You can obtain more information from your local authorities.

The batteries should be removed prior to disposal of the device and disposed of separately from the device. To protect the environment, non-rechargeable and rechargeable batteries may not be disposed of with normal household waste, but must be taken to suitable collection points. Also note the relevant statutory regulations governing the disposal of batteries.

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